

STEAMGOOGLE : INNOVATIVE MOIST HEAT TREATMENT FOR DRY EYE RELIEF



Utilizes Innovative
Heater
Technology



Safe and Effective
Dry Eye Relief



Easy Setup
and Use



Fast Warm
Up Time



Pinpointed
Moist Heat
Therapy


Relief From Symptoms of MGD and Blepharitis

Read/Watch TV During Treatment

Use Anywhere With The Included Portable Power Bank

EFFECTIVE RESULTS

The patented SteamGoggle is an easy to use at home dry eye treatment device that uses moist heat to alleviate symptoms caused by MGD and Blepharitis. The device is easy to setup (simply add water to each eyecup and press power button) and clean (place entire goggle under faucet to rinse), while also allowing the user to read/watch TV during treatment. Additionally, innovative PTC heater technology with embedded safety features allows the device to heat up fast and precisely, which allows for a high concentration of steam for each treatment. This sustained high concentration of moist heat (may take some users several treatment sessions for their eyes to acclimate to) allows the device to provide a more effective dry eye treatment compared to other methods. Use the device anywhere, anytime when using the included portable power bank with over 10 treatments per charge.

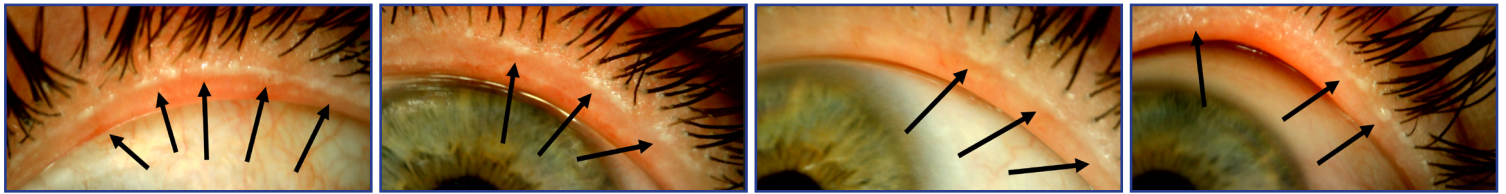

Thickened, white oil gland contents are flushed out and partially deposited on lid margin during treatment. Images were each taken following SteamGoggle treatment without any physical gland expression.

IMPORTANT INFO

- ⌚ Treatment must be performed while sitting up with face looking down or straight ahead. Do not look up or lay head down during treatment to ensure heated water is retained in eye cups.
- ⌚ Note that while looking downwards, you will experience a warmer treatment temperature versus looking straight ahead.
- ⌚ If the temperature is uncomfortable, loosen the goggles or try the treatment without adding water to experience a cooler temperature treatment. It may take several treatments for the eyes to acclimate to the moist heat.
- ⌚ Keeping the eyes open may increase the effectiveness of the treatment as blinking (consider trying some forced blinks) may help break up and remove solidified oil gland contents.
- ⌚ Use the device 1-2 times per day for 2-3 weeks, to evaluate effectiveness. Afterwards, use as needed to keep symptoms at bay. Always consult your doctor or other qualified healthcare provider for additional medical advice regarding your condition.

ABOUT

Starting in late high school and continuing through college and my working life, I have suffered through a continuous bout of symptoms related to dry eye, involving a nonstop gritty, scratchy, foreign body sensation in my left eye.

After years of finding no relief from current products/treatments on the market, (being a mechanical engineer) I created a new device that finally ended up alleviating my symptoms. The SteamGoggle differs from other products on the market, as it uses a pinpointed approach to deliver effective amounts of sustained moist heat directly to the eyes, all in a safe, easy to use at-home package.



*Power Cable Not Shown



Instructions for Use

- 1 With clean hands, begin by adding water (clean tap water is fine) to the plastic lattice inside of each eye cup by using the squirt bottle or by pouring water into the eyecup.




- 2 Tilt the goggles upside down and lightly tap with finger to remove any excess water. Some water should remain in the plastic grooves.



- 3 Connect the device to electric power, via wall plug or a portable power bank using its USB-C 12 volt output port. There should be an audible beep from power switch. Next, activate the heater by pressing the power button.

Wall Plug

Power Bank

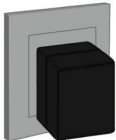

or



- 4 Allow the device to pre-heat for 1 minute 30 seconds to allow device to come to temperature. Power switch will emit an audible double-beep to notify user that pre-heating has finished.

Preheat



- 5 Again, tilt the goggles upside down and lightly tap with finger to remove any excess warmed water. Some water should remain in the plastic grooves.




- 6 Put the goggles loosely over the eyes to begin treatment. The lens may initially be cloudy, but the bottom section should shortly defog to the point where the user can look at their phone/TV.



- 7 At the end of the 8 minute 30 second treatment time, the power switch will triple-beep and shut off power. If desired, treatment can be extended by restarting power button and refilling plastic grooves with water if it has evaporated.



- 8 Clean the eyelids and device with clean water and a mild cleanser (i.e. we recommend baby soap) after each treatment session. Note that the entire goggle is water resistant (front and back), and can be placed under the faucet to easily wash away any soap remnants. Just ensure that water does not get on the power switch or power bank.

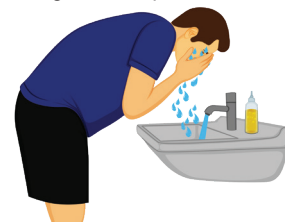

IMPORTANT: Treatment must be performed while sitting up with face looking down or straight ahead. Do not look up or lay head down during treatment to ensure heated water is retained in eyecup.